

Why **newFREEWAY**?


The Possum Freeway, trusted by thousands of satisfied users, has been dramatically re-designed to further enrich the lives of people with visual or cognitive impairment.

Using the latest technology, **newFreeway** offers enhanced levels of control and accessibility while maintaining the adaptability and reliability you can expect from a Possum product.

Environmental control enables people with physical and other impairments to control their living space.

newFREEWAY

- ✓ easy to see
- ✓ easy to hear
- ✓ **easy to use**

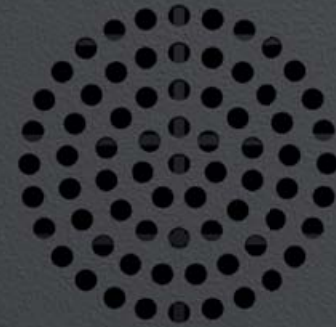

possum 

Possum Ltd is the market leader in the design, development, manufacture and installation of Electronic Assistive Technology (EAT) products and services; providing practical solutions which enhance the independence and quality of life for people with disabilities. With over 50 years' experience in Electronic Assistive Technology, we understand and respond to the needs of our clients and provide advice on solutions for a variety of differing abilities.

Possum Ltd, 8 Farmbrough Close,
Aylesbury, Bucks HP20 1DQ, UK

01296 461000 newfreeway@possum.co.uk

www.newfreeway.co.uk

possum 

NEW FREEWAY
liveyourway

Accessible Environmental Control

FREEWAY

What does newFREEWAY do?

newFreeway is an infrared (IR) and radio frequency (RF) environmental controller that provides simple remote control of a wide range of domestic appliances and entertainment devices.

The intuitive interface can be accessed in 3 easy ways, including: **Integral touch switch** | **Wired switch** | **IR Remote switch** providing a wider range of users access to environmental control and greater independence.

An integral Possum Telecare RF transmitter enables discreet activation of numerous wireless emergency communication devices including: the Possum Info Pager, Possum Neo Care phone or an existing nurse call system.

The stylish, yet robust newFreeway can be used as a standalone device or mounted on a wheelchair for operation in any room.

A large, extremely clear, backlit scanning display provides a choice of icon and / or text and can include voice announcement to further enhance the scanning functionality. A fully configurable overlay gives the user complete control over the choice of command format. newFreeway supports 3 to 32 functions that can be organised into groups for user convenience. Each function is able to generate multiple actions in sequence and can be tailored to specific room or individual user requirements.

The ability to play back specific pre-recorded spoken phrases also provides supplementary communication assistance. In addition, users can send carers pre-defined messages direct to their Info Pager.

Approvals

Product conforms to: EN60601-1 (Safety) EN60601-1-2 (EMC) following the provisions of the Medical Devices Directive.

What is different?

A technically enhanced range of features has been designed and developed to benefit improved accessibility, communication and control.

newFreeway features:

- 33% larger, brighter, backlit screen for clearer viewing
- Auto dim available to suit ambient light level
- Enhanced user control
 - Easy to operate touch switch with wipe clean surface
 - Single or two switch scanning
 - Momentary, held, manual input modes
- Scanning options
 - Linear, group, turbo
- Lighter, slim line, wipe clean casing
- 7500 mAh Li-ion battery (tested to EN62133)
- Convenient on board adjustability of:
 - Scan speed
 - Volume
 - Switch sensitivity
 - Display brightness
 - On / Off
- Easy to use PC software for programming with overlay design wizard
- Auditory scan prompt
 - Recorded Speech
 - Multi-tonal sounds
- Auditory selection confirmation
 - Recorded Speech
 - Multi-tonal sounds



Touch switch



Wired switch

Remote switch

HC8002

